



Grasp Bag Filter

BG Type

Major Applications

Highly concentrated extracts

Process water

Quality standards

- Adaptable Food Sanitation Standard
• FDA 21CFR
- Traceability by lot number

Features

- Bag construction with polypropylene filter media
- Variety of micron ratings
- Surface-treated filter material

Advantages

- High dust holding capacity, easy replacement
- Optimal micron ratings selection according to the purpose
- Reduction of fiber release

Specifications



Materials	Grades	01	05	10	25	50	75	1H	
	Micron Ratings (μ m)	1.0	5.0	10	25	50	75	100	
	Size: E.F.A. (m²)	D1, S1 : 0.25		D2, S2 : 0.5		D3, S3 : 0.09		D4, S4 : 0.15	D5, S5 : 0.25
	Size: Capacity (L)	D1, S1 : 8.0		D2, S2 : 17.5		D3, S3 : 1.5		D4, S4 : 2.5	D5, S5 : 8.0
	Media	Polypropylene							
Sealing	Polypropylene／SUS304								
Maximum Operating Temp		80 °C (176°F)							
Recommended exchange differential pressure		0.10MPa (15psi)							

*If you need further information on specifications (length, end cap type, etc.), please contact us.

Differential Pressure vs Flow Rate

Fluid: Refined Water 20~30°C (68~86°F) / Micron Ratings: 1.0~5.0 μm

Fluid: Refined Water 20~30°C (68~86°F) / Micron Ratings: 10~100 μm

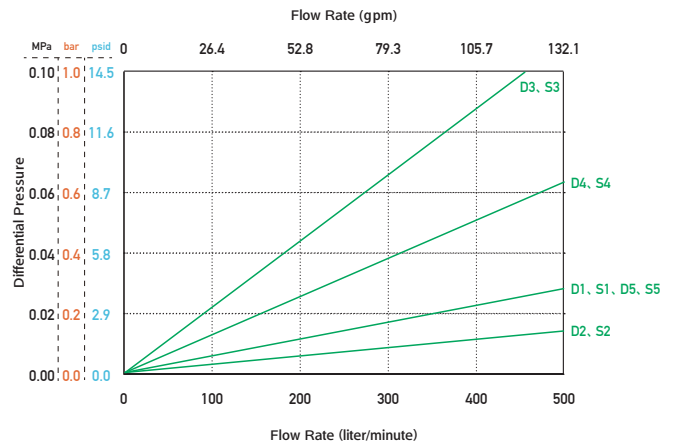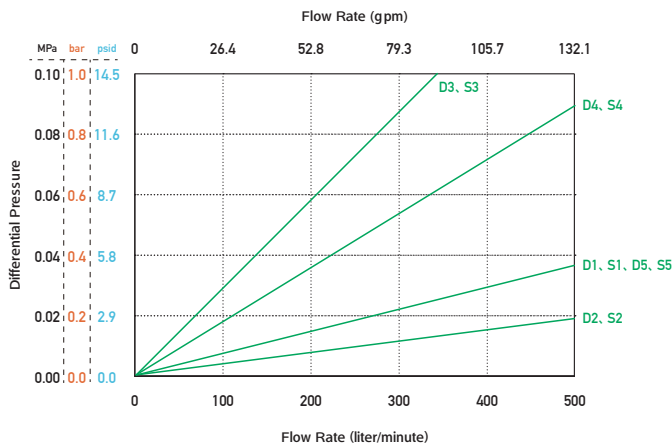

Ordering Information

Product Type	Micron Rating	Media	Sealing	Size	Packaging Code
BG	0 1	SC-	P	- D	1
	▼	▼	▼	▼	▼
	01 = 1.0 μm 05 = 5.0 μm 10 = 10 μm 25 = 25 μm 50 = 50 μm 75 = 75 μm 1H = 100 μm	P = Polypropylene	P = Polypropylene T = SUS304	Shape*1 D = Molding parts S = Snap ring	Diameter/Lenth*2 1 = 180mm / 420mm 2 = 180mm / 820mm 3 = 100mm / 230mm 4 = 100mm / 380mm 5 = 150mm / 510mm
					C = 10pcs G = 50pcs

*1 If the sealant is SUS304, please select snap ring shape.

*2 If the sealant is polypropylene with snap ring shape, there is no setting for diameter/overall length 3, 4, and 5.

*The contents of the catalog are subject to change without notice.

*The performance data listed in the catalog are Typical values obtained under specific conditions based on our tests.

ROKITECHNO MIRAI CO., LTD.

6-20-12, Minami-Oi, Shinagawa-ku Tokyo, 140-0013 Japan

TEL: +81-3-5764-1131 FAX: +81-3-5764-0681

www.rokitechno.com

For our technical information, please click here.▼



Manufacturing is based on our Quality Management Systems that meet ISO9001 standards.

Scope

Design, Development, manufacture, and sales of filter cartridges, housings and filtration equipment.



4th Issue
BG210301LE